FTGN-2205

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ENGLISH

Read these instructions for use carefully and keep them for future reference. Please keep them for later use, make them available to other users and follow the instructions.

WARNING

- The appliance is intended for domestic/private use only, not for commercial use.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children must not play with the appliance.
- Cleaning and user maintenance must not be carried out by children without supervision.
- If the mains connection cable of this appliance is damaged, it must be disposed of. If it cannot be removed, the appliance must be disposed of.
- The appliance has a hot surface. People who are not sensitive to heat must be careful when using the appliance.
- Do not insert any objects into the appliance openings and do not insert anything into the moving parts. Make sure that the moving parts can always move freely.
- Do not open or repair the appliance under any circumstances, as otherwise it will no longer function properly. Failure to do so will invalidate the warranty.

Dear customer,

We are delighted that you have chosen a product from our range. Our name stands for high-quality and thoroughly tested fitness and wellness products.

With kind recommendation,

Your Fitagon team

Contents

1. scope of delivery.....	12	7. commissioning.....	16
2. explanation of symbols.....	12	8. application.....	16
3. Get to know.....	13	9. cleaning and storage.....	17
4. intended use.....	13	10. disposal.....	17
5. warnings and safety instructions.....	14	11. warranty and service.....	18
6. device description.....	15		

1. Scope of delivery

Check the device to ensure that the packaging is intact and that the contents are complete. Before use, ensure that the appliance and accessories do not show any visible signs of damage and all packaging material is removed. If in doubt, do not use it and contact to your dealer or to the specified service address.

1 x head and body massager


1 x charging station

1 x Charging cable

1 x These instructions for use

2. Explanation of symbols

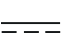The following symbols are used in the operating instructions and on the rating plate:

	WARNING Warning of risk of injury or danger to your health.
	ATTENTION Safety note on possible damage to the device/accessories.
	Note Reference to important information.
	Safety isolating transformer, short-circuit proof
	Use only in closed rooms.
	Dispose of packaging in an environmentally friendly manner
	Disposal in accordance with the Waste Electrical and Electronic Equipment (WEEE) Directive
	Manufacturer
	CE marking This product meets the requirements of the applicable European and national directives.
	Read instructions for use
	Direct current Device is only suitable for direct current
	Energy efficiency level 6
IP20	Protected against solid foreign objects, 12.5 mm diameter and larger
	The products demonstrably meet the requirements of the technical regulations of the EAEU
	Polarity of the DC power connection

3. Get to know the device

The Fitagon Head and Body Massager is a waterproof and cordless companion that relieves headaches and muscle tension while promoting hair growth. Particularly suitable for those who suffer from neck pain, this massager can be safely used while showering to relieve itching and improve circulation.

Its versatile application allows for a full body massage and even provides gentle body exfoliation. The 4 massage heads are equipped with a total of 27 soft silicone massage points which, unlike conventional massage claws, do not get tangled in the hair.

Easily portable and ideal for everyday use at home or on the go, the massage heads can be easily removed and cleaned. With its elegant design, the Fitagon head and body massager is an ideal gift choice for friends and family.

The head massager is for personal use only, not for medical or commercial use.

4. Intended use

The Fitagon head and body massager relieves pain and promotes relaxation. It is not a substitute for medical treatment. Do not use the massager if one or more of the following warnings apply to you. If you are unsure whether the massager is suitable for you, consult your doctor.

WARNING

Use the massager:

- not with animals,
- not in the case of a pathological change or injury in the body region to be massaged (e.g. herniated disc, open wound),
- Never use it for massages in the heart area if you have a pacemaker and consult your doctor first for other areas of the body,
- Never use on swollen, burnt, inflamed or injured areas of skin or body,
- not for lacerations, capillaries, varicose veins, acne, couperose veins, herpes or other skin diseases,
- not on the head and hard parts of the body (e.g. bones and joints),
- not during pregnancy,
- not on the face (eyes), larynx or other particularly sensitive parts of the body,
- never while you are sleeping,
- not after taking medication or alcohol (limited ability to perceive!),
- not on creamed skin (e.g. warming ointment),
- not in the motor vehicle.

Consult your doctor before using the device, in particular:

- if you are suffering from a serious illness or have had an operation,
- if you wear a pacemaker, implants or other aids,
- for thromboses,
- for diabetes,
- for pain of unknown cause.

The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage caused by improper or careless use.

5. Warnings and safety instructions

Read these instructions for use carefully! Failure to observe the following instructions may result in personal injury or damage to property. Keep these instructions for use in a safe place and make them available to other users. Hand over these instructions for use when passing on the appliance.

WARNING

Keep children away from packaging material. There is a risk of suffocation.

Electric shock

WARNING

As with any electrical appliance, this appliance must be used with care and caution to avoid the risk of electric shock.

Therefore, operate the massager:

- only on the mains voltage specified on the mains adapter,
- never if the massager or the accessories are visibly damaged,
- not during a thunderstorm.

In the event of defects or malfunctions, switch the massager off immediately and disconnect the massager from the power supply. Do not pull on the mains cable or on the massager to unplug the mains adapter from the socket. Hold or carry

Never connect the appliance to the mains cable. Keep a distance between the cables and hot surfaces.

Do not pinch, kink or twist the mains cable. Do not use needles or sharp objects pierce it.

Ensure that the massager, the remote control, the mains adapter and the mains cable do not come into contact with water, steam or other liquids. Therefore, only use the massager in dry indoor areas (e.g. never in the bathtub or sauna).

Do not use the device if there is visible damage to the device or accessories. Do not subject the device to impacts and do not drop it.

Repair

WARNING

- Repairs to electrical appliances may only be carried out by qualified personnel. Due to improper

Repairs can cause considerable danger to the user. Contact customer service or an authorized dealer for repairs.

- Do not misuse the mains cable.

Fire hazard

WARNING

In the event of improper use or non-compliance with these instructions for use there may be a risk of fire!

Therefore, operate the massager

- never under a cover, such as a blanket, cushion, ...
- Never use near gasoline or other highly flammable substances.

Handling

WARNING


The appliance must be switched off and unplugged after each use and before cleaning.

- Do not sit, lie or stand with your full weight on the moving parts of the appliance and do not place any weight on them.


Place objects on the appliance.

- Protect the device from high temperatures.

6. Device description

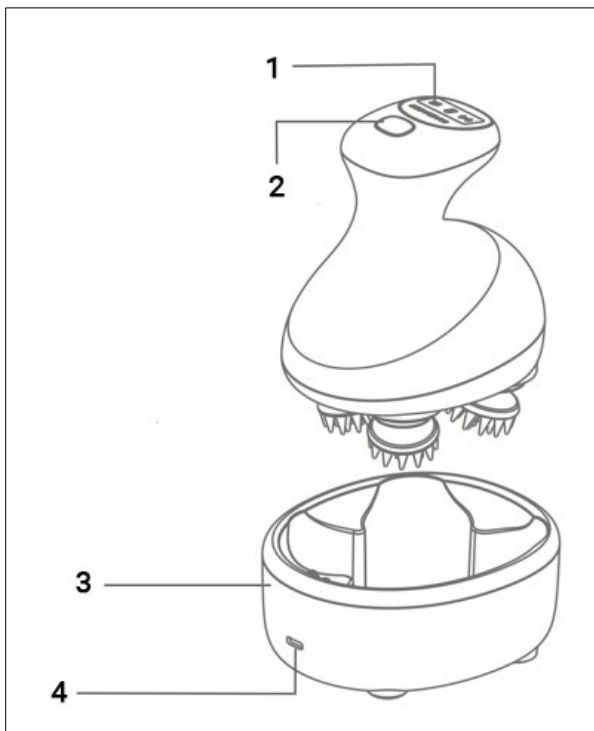1. display symbols
2. switch button
3. charging station
4. type-C charging input

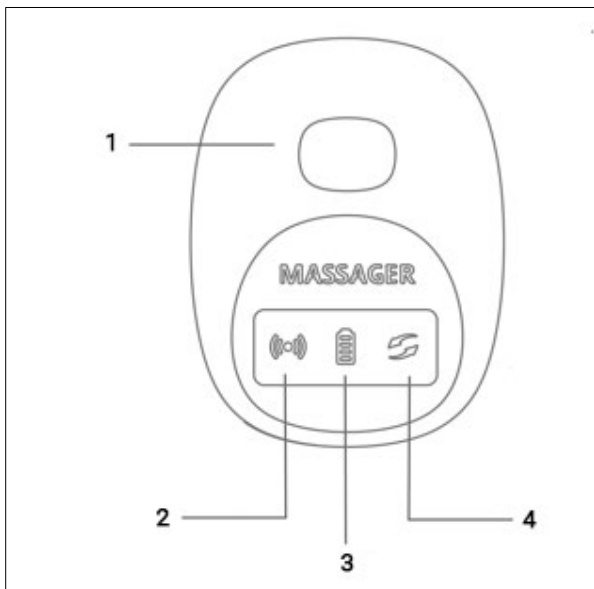1. switch button
2. vibration
3. power
4. massage

7. Commissioning

- Remove the packaging.
- Check the device, charging station and charging cable for damage.
- Fully charge the device before using it for the first time. Please use a 5V/1A charger for charging.
- It takes about 2-3 hours to fully charge, but due to the ambient temperature, the charging time and the usage time of the product may vary.

8. Application

ON/OFF

Pressing and holding the switch button switches the device on/off. Pressing the switch button briefly switches the operating mode.

Modes

Mode 1: Slow speed with red light.

Mode 2: Fast speed with red light.

Mode 3: Mixing slow and fast speed with red light.

Mode 4: Slow speed without red light.

Mode 5: Fast speed without red light.

Mode 6: Mixing slow and fast speed without red light.

Light indicator

When the device is switched on, the indicator light flashes white and the red-light flashes.

The red-light flashes quickly when the battery is low, as well as during the charging process.

The white light is always on when the battery is fully charged. When the speed is slow,

the indicator flashes slowly, at high speeds the indicator light flashes quickly. The indicator light changes with the speed.

The white light of the massage icon flashes slowly when the speed is slow, the white light of the massage icon flashes quickly when the speed is high, and the white light changes with the speed when the slow and fast speeds are mixed.

Double-click the vibration button to activate the vibration function. At this point, the white light of the vibration icon flashes three times. Press the button twice again to switch off the vibration function.

Battery indicator light

The battery indicator lights up white when the device is switched on, the red light flashes quickly when the battery is low and the red light flashes during the charging process. The white light is always on when the battery is fully charged.

Note

- The massage should be perceived as pleasant and relaxing at all times. Stop the massage if you find it painful or unpleasant.

- Do not use the device for longer than 15-30 minutes. If the massage lasts longer, overstimulation of the muscles can lead to tension instead of relaxation.
- The device has an automatic switch-off function that is set to a maximum operating time of 10 minutes.

Notes on handling rechargeable batteries

- If liquid from a battery cell comes into contact with skin or eyes, rinse the affected area with water and seek medical attention.
- Danger of swallowing! Small children could swallow batteries and choke on them. Therefore, keep batteries out of reach of small children!
- Protect batteries from excessive heat.
- Explosion hazard! Do not throw batteries into a fire.
- Do not disassemble, open or crush batteries.
- Only use chargers listed in the instructions for use.
- Batteries must be charged correctly before use. Always follow the manufacturer's instructions or the information in these instructions for use for correct charging.
- Charge the battery fully before using it for the first time.
- To achieve the longest possible battery life, fully charge the battery at least twice a year.

9. Cleaning and storage

WARNING

- Always disconnect the appliance from the power supply before cleaning.

Cleaning

- Only clean the appliance in the manner specified. Never allow any liquid to enter the appliance or accessories.
- Smaller stains can be removed with a cloth and possibly a little liquid mild detergent. Do not use cleaning agents containing solvents.
- Do not use the appliance again until it is completely dry.

Storage

If you do not use the device for a longer period of time, we recommend that you store it in its original packaging in a dry environment and without weight.

10. Disposal

ATTENTION

In the interests of environmental protection, the device must not be disposed of with household waste at the end of its service life. It can be disposed of at the appropriate collection points in your country. Dispose of the device in accordance with the Waste Electrical and Electronic Equipment (WEEE) Directive. If you have any questions, please contact the local authority responsible for disposal.

Battery disposal

You must dispose of used, fully discharged batteries via specially marked collection collection containers, special waste collection points or via the electrical retailer. You are legally obliged to dispose of the batteries. To remove the battery from the device remove the battery from the device, pry open the cover above the light ring with a screwdriver with a screwdriver. Remove the battery from the inside of the housing. You will find these symbols on batteries containing harmful substances:

Pb = battery contains lead,

Cd = battery contains cadmium,

Hg = battery contains mercury.

11. Warranty and service

Fitagon GmbH, Technoparkstrasse 1, CH-8406 Winterthur grants a warranty for this product under the following conditions and to the extent described below.

The following warranty conditions do not affect the seller's statutory warranty obligations arising from the purchase contract with the buyer. The guarantee also applies without prejudice to mandatory statutory liability provisions.

Fitagon guarantees the faultless functionality and completeness of this product. The warranty period is 2 years from the date of purchase of the new, unused product by the buyer.

This warranty only applies to products purchased by the buyer as a consumer and used exclusively for personal purposes within the scope of domestic use. Swiss law applies.

If, during the warranty period, this product proves to be incomplete or defective in function in accordance with the following provisions, Fitagon will provide a replacement delivery or repair free of charge in accordance with these warranty conditions.

If the buyer wishes to report a warranty claim, they should first contact Fitagon Customer Service:

Fitagon GmbH

Mail: info@fitagon.ch

Tel: +41 52 202 57 70

The buyer will then receive more detailed information on how to handle the warranty claim, e.g. where to send the product free of charge and which documents are required.

A claim under the guarantee can only be considered if the buyer

- a copy of the invoice/purchase receipt and
- the original product

Fitagon or an authorized Fitagon partner.

Expressly excluded from this guarantee are

- Wear and tear resulting from normal use or consumption of the product;
- Accessories supplied with this product that wear out or become damaged during normal use. are consumed (e.g. batteries, rechargeable batteries, cuffs, seals, electrodes, illuminants, attachments, inhaler accessories);
- Products that are used improperly and/or contrary to the provisions of the operating instructions have been used, cleaned, stored or maintained, as well as products that have been used, cleaned, stored or maintained by the buyer or not by Fitagon have been opened, repaired or converted;
- Damage that occurs during transport between the manufacturer and the customer or between the service center and the customer;
- Products that were purchased as 2nd choice items or as used items;

- consequential damage resulting from a defect in this product (however, in this case Claims arising from product liability or other mandatory statutory liability provisions exist).
Repairs or a complete replacement do not extend the warranty period under any circumstances.

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