# FTGN-2207

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# ENGLISH

Read these instructions for use carefully and keep them for future reference. Please keep them for later use, make them available to other users and follow the instructions.

## WARNING

- The appliance is intended for domestic/private use only, not for commercial use.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children must not play with the appliance.
- Cleaning and user maintenance must not be carried out by children without supervision.
- If the mains connection cable of this appliance is damaged, it must be disposed of. If it cannot be removed, the appliance must be disposed of.
- The appliance has a hot surface. People who are not sensitive to heat must be careful when using the appliance.
- Do not insert any objects into the appliance openings and do not insert anything into the moving parts. Make sure that the moving parts can always move freely.
- Do not open or repair the appliance under any circumstances, as otherwise it will no longer function properly. Failure to do so will invalidate the warranty.

### Dear customer,

We are delighted that you have chosen a product from our range. Our name stands for high-quality and thoroughly tested fitness and wellness products.

With kind recommendation,

Your Fitagon team

## Contents

1. scope of delivery.....	12	7. commissioning.....	16
2. explanation of symbols.....	12	8. application.....	16
3. Get to know.....	13	9. cleaning and storage.....	17
4. intended use.....	13	10. disposal.....	17
5. warnings and safety instructions.....	14	11. warranty and service.....	18
6. device description.....	15		

# 1. Scope of delivery

Check the device to ensure that the packaging is intact and that the contents are complete. Before use, ensure that the appliance and accessories do not show any visible signs of damage and all packaging material is removed. If in doubt, do not use it and contact to your dealer or to the specified service address.

1 x neck stretcher with heat and impulse massage with remote control


1 x electrical impulse patch (with connection cable and two patches)

1 x mains adapter

1 x These instructions for use

# 2. Explanation of symbols

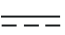The following symbols are used in the operating instructions and on the rating plate:

	<b>WARNING</b> Warning of risk of injury or danger to your health.
	<b>ATTENTION</b> Safety note on possible damage to the device/accessories.
	<b>Note</b> Reference to important information.
	Safety isolating transformer, short-circuit proof
	Use only in closed rooms.
	Dispose of packaging in an environmentally friendly manner
	Disposal in accordance with the Waste Electrical and Electronic Equipment (WEEE) Directive
	Manufacturer
	<b>CE marking</b> This product meets the requirements of the applicable European and national directives.
	Read instructions for use
	<b>Direct current</b> Device is only suitable for direct current
	Energy efficiency level 6
<b>IP20</b>	Protected against solid foreign objects, 12.5 mm diameter and larger
	The products demonstrably meet the requirements of the technical regulations of the EAEU
	Polarity of the DC power connection

### 3. Get to know the device

The Fitagon neck stretcher offers an innovative combination of heat, pulse massage and air traction to promote neck health. With three automatic modes and dynamic traction that can lift the neck by up to 10 cm, it offers a variety of benefits. The heat treatment relieves muscle stiffness and promotes circulation, while the electrical stimulation penetrates deep into the skin to soothe muscles. This device also provides a full body massage thanks to the low frequency impulse electrotherapy and is an effective way to relieve neck discomfort.

The neck stretcher device is for personal use only, not for medical or commercial use.

### 4. Intended use

This device mobilizes and stretches the cervical spine while massaging the neck muscles. It is not a substitute for medical treatment. Do not use the massager if one or more of the following warnings apply to you. If you are unsure whether the massager is suitable for you, consult your doctor.

#### WARNING

Use the massager:

- not with animals,
- not in the case of a pathological change or injury in the body region to be massaged (e.g. herniated disc, open wound),
- Never use it for massages in the heart area if you have a pacemaker and consult your doctor first for other areas of the body,
- Never use on swollen, burnt, inflamed or injured areas of skin or body,
- not for lacerations, capillaries, varicose veins, acne, couperose veins, herpes or other skin diseases,
- not on the head and hard parts of the body (e.g. bones and joints),
- not during pregnancy,
- not on the face (eyes), larynx or other particularly sensitive parts of the body,
- never while you are sleeping,
- not after taking medication or alcohol (limited ability to perceive!),
- not on creamed skin (e.g. warming ointment),
- not in the motor vehicle.

Consult your doctor before using the device, in particular:

- if you are suffering from a serious illness or have had an operation,
- if you wear a pacemaker, implants or other aids,
- for thromboses,
- for diabetes,
- for pain of unknown cause.

The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage caused by improper or careless use.

## 5. Warnings and safety instructions

Read these instructions for use carefully! Failure to observe the following instructions may result in personal injury or damage to property. Keep these instructions for use in a safe place and make them available to other users. Hand over these instructions for use when passing on the appliance.

### WARNING

Keep children away from packaging material. There is a risk of suffocation.

Electric shock

### WARNING

As with any electrical appliance, this appliance must be used with care and caution to avoid the risk of electric shock.

Therefore, operate the massager:

- only on the mains voltage specified on the mains adapter,
- never if the massager or the accessories are visibly damaged,
- not during a thunderstorm.

In the event of defects or malfunctions, switch the massager off immediately and disconnect the massager from the power supply. Do not pull on the mains cable or on the massager to unplug the mains adapter from the socket. Hold or carry

Never connect the appliance to the mains cable. Keep a distance between the cables and hot surfaces.

Do not pinch, kink or twist the mains cable. Do not use needles or sharp objects pierce it.

Ensure that the massager, the remote control, the mains adapter and the mains cable do not come into contact with water, steam or other liquids. Therefore, only use the massager in dry indoor areas (e.g. never in the bathtub or sauna).

Do not use the device if there is visible damage to the device or accessories. Do not subject the device to impacts and do not drop it.

Repair

### WARNING

- Repairs to electrical appliances may only be carried out by qualified personnel. Due to improper

Repairs can cause considerable danger to the user. Contact customer service or an authorized dealer for repairs.

- Do not misuse the mains cable.

Fire hazard

### WARNING

In the event of improper use or non-compliance with these instructions for use there may be a risk of fire!

Therefore, operate the massager

- never under a cover, such as a blanket, cushion, ...
- Never use near gasoline or other highly flammable substances.

#### Handling

#### WARNING


The appliance must be switched off and unplugged after each use and before cleaning.

- Do not sit, lie or stand with your full weight on the moving parts of the appliance and do not place any weight on them.

Place objects on the appliance.

- Protect the device from high temperatures.

## 6. Device description

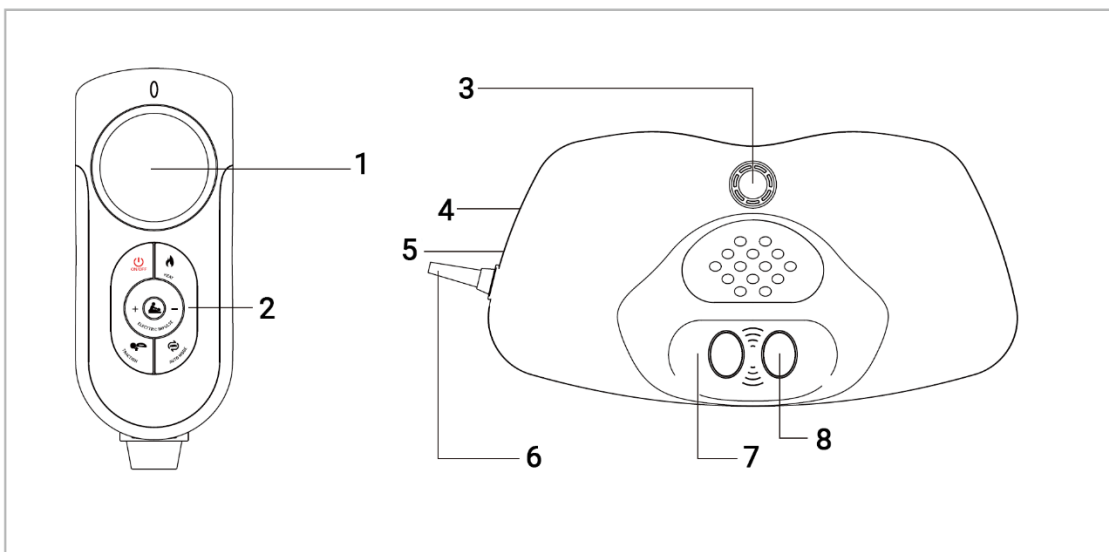1. screen
2. operating buttons
3. power light
4. mains connection
5. output for electrical pulses
6. control cable
7. heating zone
8. electrical, pulse-generating foil

## 7. Commissioning

- Remove the packaging.
- Check the device and the mains adapter for damage.
- Insert the adapter into the device and connect it to the power supply. The remote control always stays on when the device is charged
- Lay the mains cable so that it cannot be tripped over.

# 8. Application

## ON/OFF

Press the On/Off button to switch the appliance on or off. The operating time on the display is 15 minutes.

## Pressure function

Press the  button and the pressure rises and opens. The pressure is divided into three modes: A/B/AB

A - Comfort mode

B - Boost tension

AB - Intelligent movement mode

Press the Pulse button to change the mode.

## Electrical pulse

Press the  button to switch on the electrical pulse. There are three modes: A/B/AB, including the plus and minus buttons to increase or decrease the intensity.

A - Pulse

B - Acupuncture

AB - Intelligent pulse

There are 12 different levels for each mode

When using the device, a curve on the screen refers to a level of intensity of the electrical impulse. If the intensity is increased from 7 to 12, the screen will still show curves 1 to 6, but flashing.

Press the Pulse button to change the mode.

## Heat

Press the  button to switch on the heating mode. The heating mode is divided into three intensities: A/B/AB.

A - low temperature

B - medium temperature

AB - high temperature


Press the heating button to switch between the temperatures.

## Auto mode

Press the  button to activate automatic mode. In automatic mode, pressure, electrical impulse, heat and acupuncture are combined. Press the auto button to switch between the temperatures.

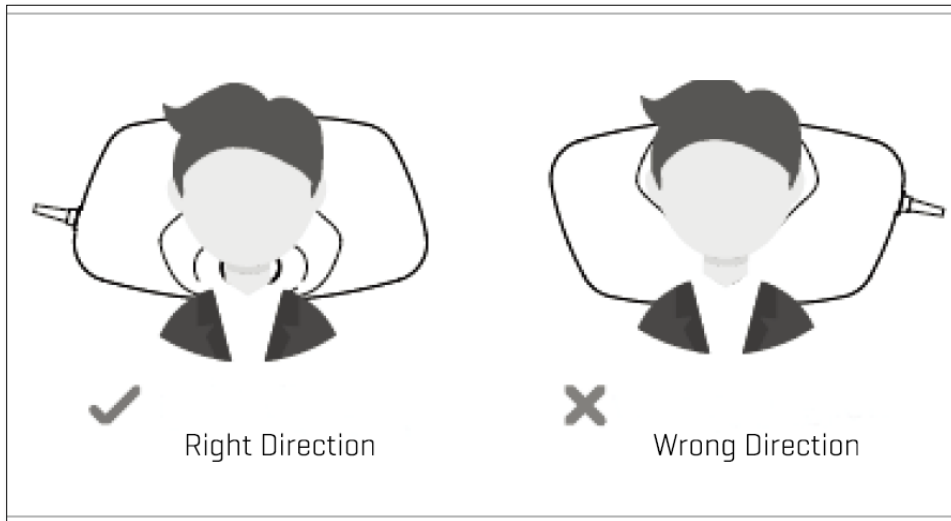Place the device directly under the neck, plug it into the power supply and use can begin.

### Note

- The massage should be pleasant and relaxing at all times. Stop the massage or change the pressure using the handles if you find the massage painful or unpleasant.
- Do not use the device for longer than 15-30 minutes. If the massage lasts longer, overstimulation of the muscles can lead to tension instead of relaxation
- The appliance has an automatic switch-off function that is set to a maximum operating time of 15 minutes

## 9. Cleaning and storage

### WARNING

- Always disconnect the appliance from the power supply before cleaning.

### Cleaning

- Only clean the appliance in the manner specified. Never allow any liquid to enter the appliance or accessories.
- Smaller stains can be removed with a cloth and possibly a little liquid mild detergent. Do not use cleaning agents containing solvents.
- Do not use the appliance again until it is completely dry.

### Storage

If you do not use the device for a longer period of time, we recommend that you store it in its original packaging in a dry environment and without weight.

## 10. Disposal

### ATTENTION

In the interests of environmental protection, the device must not be disposed of with household waste at the end of its service life. It can be disposed of at the appropriate collection points in your country. Dispose of the device in accordance with the Waste Electrical and Electronic Equipment (WEEE) Directive. If you have any questions, please contact the local authority responsible for disposal.

# 11. Warranty and service

Fitagon GmbH, Technoparkstrasse 1, CH-8406 Winterthur grants a warranty for this product under the following conditions and to the extent described below.

**The following warranty conditions do not affect the seller's statutory warranty obligations arising from the purchase contract with the buyer. The guarantee also applies without prejudice to mandatory statutory liability provisions.**

Fitagon guarantees the faultless functionality and completeness of this product. The warranty period is 2 years from the date of purchase of the new, unused product by the buyer.

This warranty only applies to products purchased by the buyer as a consumer and used exclusively for personal purposes within the scope of domestic use. Swiss law applies.

If, during the warranty period, this product proves to be incomplete or defective in function in accordance with the following provisions, Fitagon will provide a replacement delivery or repair free of charge in accordance with these warranty conditions.

**If the buyer wishes to report a warranty claim, they should first contact Fitagon Customer Service:**

**Fitagon GmbH**

**Mail: [info@fitagon.ch](mailto:info@fitagon.ch)**

**Tel: +41 52 202 57 70**

The buyer will then receive more detailed information on how to handle the warranty claim, e.g. where to send the product free of charge and which documents are required.

A claim under the guarantee can only be considered if the buyer

- a copy of the invoice/purchase receipt and
- the original product

Fitagon or an authorized Fitagon partner.

Expressly excluded from this guarantee are

- Wear and tear resulting from normal use or consumption of the product;
- Accessories supplied with this product that wear out or become damaged during normal use. are consumed (e.g. batteries, rechargeable batteries, cuffs, seals, electrodes, illuminants, attachments, inhaler accessories);
- Products that are used improperly and/or contrary to the provisions of the operating instructions have been used, cleaned, stored or maintained, as well as products that have been used, cleaned, stored or maintained by the buyer or not by Fitagon have been opened, repaired or converted;
- Damage that occurs during transport between the manufacturer and the customer or between the service center and the customer;
- Products that were purchased as 2nd choice items or as used items;
- consequential damage resulting from a defect in this product (however, in this case Claims arising from product liability or other mandatory statutory liability provisions exist).

Repairs or a complete replacement do not extend the warranty period under any circumstances.

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