FTGN-2201

FTGN-2202

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ENGLISH

Read these instructions for use carefully and keep them for future reference. Please keep them for later use, make them available to other users and follow the instructions.

WARNING

- The appliance is intended for domestic/private use only, not for commercial use.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children must not play with the appliance.
- Cleaning and user maintenance must not be carried out by children without supervision.
- If the mains connection cable of this appliance is damaged, it must be disposed of. If it cannot be removed, the appliance must be disposed of.
- The appliance has a hot surface. People who are not sensitive to heat must be careful when using the appliance.
- Do not insert any objects into the appliance openings and do not insert anything into the moving parts. Make sure that the moving parts can always move freely.
- Do not open or repair the appliance under any circumstances, as otherwise it will no longer function properly. Failure to do so will invalidate the warranty.

Dear customer,

We are delighted that you have chosen a product from our range. Our name stands for high-quality and thoroughly tested fitness and wellness products.

With kind recommendation,

Your Fitagon team

Contents

| | | | |
|--|----|-----------------------------------|----|
| 1. scope of delivery..... | 11 | 7. commissioning..... | 14 |
| 2. explanation of symbols..... | 11 | 8. application..... | 16 |
| 3. Get to know..... | 12 | 9. cleaning and storage..... | 17 |
| 4. intended use..... | 12 | 10. disposal..... | 17 |
| 5. warnings and safety instructions..... | 13 | 11. technical specifications..... | 17 |
| 6. device description..... | 14 | 12. warranty and service..... | 18 |

1. Scope of delivery

Check the device to ensure that the packaging is intact and that the contents are complete. Before use, ensure that the appliance and accessories do not show any visible signs of damage and all packaging material is removed. If in doubt, do not use it and contact to your dealer or to the specified service address.


1 x Shiatsu Massager Pro

1 x mains adapter

1 x This instruction manual

2. Explanation of symbols

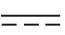The following symbols are used in the operating instructions and on the rating plate:

| | |
|---|--|
|  | WARNING Warning of risk of injury or danger to your health. |
|  | ATTENTION Safety note on possible damage to the device/accessories. |
|  | Note Reference to important information. |
|  | Safety isolating transformer, short-circuit proof |
|  | Use only in closed rooms. |
|  | Dispose of packaging in an environmentally friendly manner |
|  | Disposal in accordance with the Waste Electrical and Electronic Equipment (WEEE) Directive |
|  | Manufacturer |
|  | CE marking This product meets the requirements of the applicable European and national directives. |
|  | Read instructions for use |
|  | Direct current Device is only suitable for direct current |
|  | Energy efficiency level 6 |
| IP20 | Protected against solid foreign objects, 12.5 mm diameter and larger |
|  | The products demonstrably meet the requirements of the technical regulations of the EAEU |
|  | Polarity of the DC power connection |

3. Get to know the device

Shiatsu, which originated in Japan around a century ago, is a form of body massage based on the principles of traditional principles of traditional Chinese medicine. The basis of this technique is the system of energy channels, also known as the meridian system, in the human body. The main aim of a Shiatsu massage is to promote the physical, emotional and mental well-being of the person being treated. This is achieved by the therapist releasing energetic blockages and congestion in the energy channels and stimulating the body's self-regulating powers. To achieve this goal, the shiatsu therapist applies pressure in flowing movements along the energy channels, the meridians. Although Shiatsu is literally translated as "finger pressure", this technique fingers, this technique also uses the balls of the hands, elbows and knees. The therapist uses his or her entire body weight to generate the necessary pressure.

Your Shiatsu massage device contains rotating massage heads that replicate the pressure and kneading movements of a Shiatsu massage device. With this electric Shiatsu massager, you can give yourself or others a soothing massage effectively and without help. Massages can have a relaxing or stimulating effect and are often used to relieve muscle tension, pain and fatigue. The device provides a powerful, intensive massage for the neck, back, shoulders and thighs. The Shiatsu massager is for personal use only, not for medical or commercial use. In addition, the Shiatsu massager offers you the following advantages:

- simple operation,
- clockwise and counter-clockwise rotation,
- heat and light function,
- 3 speed settings.

4. Intended use

This device is intended for massaging the neck, back, shoulders and thighs of the human body. It is not a substitute for medical treatment. Do not use the massager if one or more of the following warnings apply to you. If you are unsure whether the massager is suitable for you, consult your doctor.

WARNING

Use the massager:

- not with animals,
- not in the case of a pathological change or injury in the body region to be massaged (e.g. herniated disc, open wound),
- Never use it for massages in the heart area if you have a pacemaker and consult your doctor first for other areas of the body,
- Never use on swollen, burnt, inflamed or injured areas of skin or body,
- not for lacerations, capillaries, varicose veins, acne, couperose veins, herpes or other skin diseases,
- not on the head and hard parts of the body (e.g. bones and joints),
- not during pregnancy,
- not on the face (eyes), larynx or other particularly sensitive parts of the body,
- never while you are sleeping,
- not after taking medication or alcohol (limited ability to perceive!),
- not on creamed skin (e.g. warming ointment),
- not in the motor vehicle.

Consult your doctor before using the device, in particular:

- if you are suffering from a serious illness or have had an operation,
- if you wear a pacemaker, implants or other aids,
- for thromboses,
- for diabetes,
- for pain of unknown cause.

The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage caused by improper or careless use.

5. Warnings and safety instructions

Read these instructions for use carefully! Failure to observe the following instructions may result in personal injury or damage to property. Keep these instructions for use in a safe place and make them available to other users. Hand over these instructions for use when passing on the appliance.

WARNING

Keep children away from packaging material. There is a risk of suffocation.

Electric shock

WARNING

As with any electrical appliance, this appliance must be used with care and caution to avoid the risk of electric shock.

Therefore, operate the massager:

- only on the mains voltage specified on the mains adapter,
- never if the massager or the accessories are visibly damaged,
- not during a thunderstorm.

In the event of defects or malfunctions, switch the massager off immediately and disconnect the massager from the power supply. Do not pull on the mains cable or on the massager to unplug the mains adapter from the socket. Hold or carry

Never connect the appliance to the mains cable. Keep a distance between the cables and hot surfaces.

Do not pinch, kink or twist the mains cable. Do not use needles or sharp objects pierce it.

Ensure that the massager, the remote control, the mains adapter and the mains cable do not come into contact with water, steam or other liquids. Therefore, only use the massager in dry indoor areas (e.g. never in the bathtub or sauna).

Do not use the device if there is visible damage to the device or accessories. Do not subject the device to impacts and do not drop it.

Repair

WARNING

- Repairs to electrical appliances may only be carried out by qualified personnel. Due to improper

Repairs can cause considerable danger to the user. Contact customer service or an authorized dealer for repairs.

- Do not misuse the mains cable.

Fire hazard

WARNING

In the event of improper use or non-compliance with these instructions for use there may be a risk of fire!

Therefore, operate the massager

- never under a cover, such as a blanket, cushion, ...
- Never use near gasoline or other highly flammable substances.

Handling

WARNING


The appliance must be switched off and unplugged after each use and before cleaning.

- Do not sit, lie or stand with your full weight on the moving parts of the appliance and do not place any weight on them.

Place objects on the appliance.

- Protect the device from high temperatures.

6. Device description

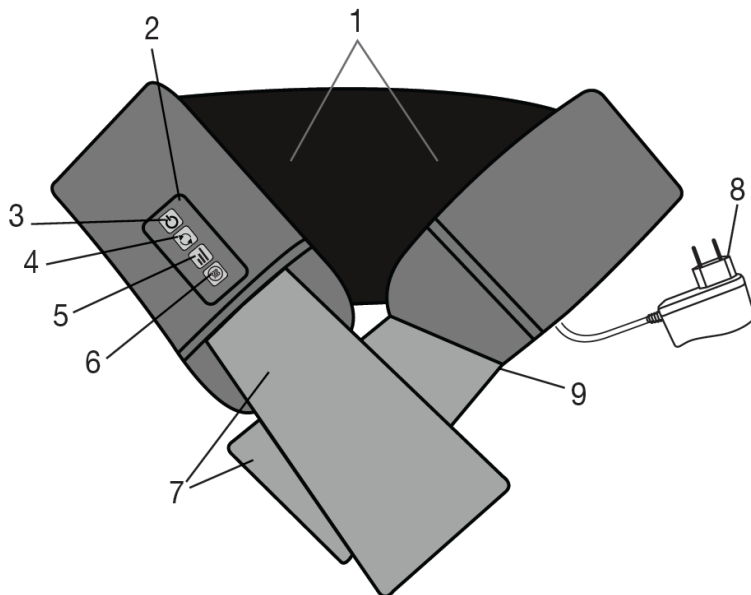1. 16x massage heads (8 heated and illuminated massage heads), rotating in pairs
2. Control element
3. -button
4. - button (anti-clockwise/clockwise rotation of the massage heads)
5. -button (3 intensity levels)
6. -button (heat function)
7. Handles
8. Mains adapter
9. Mains connection cable (cable is secured to the side)

7. Commissioning

- Remove the packaging.
- Check the device and the mains adapter for damage.
- Insert the connecting plug into the socket on the Shiatsu massage device.
- Lay the mains cable so that it cannot be tripped over.

8. Application

Switch on the Shiatsu massager with the  button. The massage heads start to rotate. You can select the massage direction by pressing the  button. If you want to switch on the heat function press the  button. If you press the  button again, the heat function switches off again. You can choose between 3 intensity levels using the  button. The device starts automatically at the lowest intensity level. To switch the Shiatsu massager off again, press the  button.

Massage variations


Neck massage:

- Place the Shiatsu massager on your neck.
- Set the contact pressure to your preferred setting using the handle.


Back massage (upper back area):

- Place the Shiatsu massager around your upper back area.
- Set the contact pressure to your preferred setting using the handle.

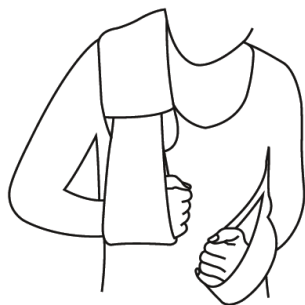Back massage (shoulder blade):

- Place the Shiatsu massager diagonally around your left or right shoulder as shown.
- Set the contact pressure to your preferred setting using the handle.

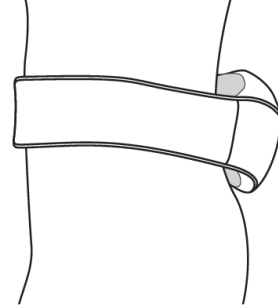Back massage (lower back area):

- Place the Shiatsu massager around your hips as shown.
- Set the contact pressure to your preferred setting using the handle.

Note

- The massage should be pleasant and relaxing at all times. Stop the massage or change the pressure using the handles if you find the massage painful or unpleasant.

- Do not use the device for longer than 15-30 minutes. If the massage lasts longer, overstimulation of the muscles can lead to tension instead of relaxation

- The appliance has an automatic switch-off function that is set to a maximum operating time of 15 minutes

9. Cleaning and storage

WARNING

- Always disconnect the appliance from the power supply before cleaning.

Cleaning

- Only clean the appliance in the manner specified. Never allow any liquid to enter the appliance or accessories.
- Smaller stains can be removed with a cloth and possibly a little liquid mild detergent. Do not use cleaning agents containing solvents.
- Do not use the appliance again until it is completely dry.

Storage

If you do not use the device for a longer period of time, we recommend that you store it in its original packaging in a dry environment and without weight.

10. Disposal

ATTENTION

In the interests of environmental protection, the device must not be disposed of with household waste at the end of its service life. It can be disposed of at the appropriate collection points in your country. Dispose of the device in accordance with the Waste Electrical and Electronic Equipment (WEEE) Directive. If you have any questions, please contact the local authority responsible for disposal.

11. Technical data

| | |
|--|---|
| Mains operation: Entrance Output | 100–240 V~, 50/60 Hz, 0,5 A 12,6 V  1,0 A, 12,6 W |
| Average efficiency in operation | ≥ 87,38 % |
| Efficiency at low load (10%) | ≥ 85,66 % |
| Power consumption at no load | ≤ 0,072 W |
| Mass | 128 x 15,7 x 14 cm |
| Weight | ca. 1,7 kg |

What to do if problems occur

| Problem | Cause | Solution |
|------------------------------|---------------------------------------|---|
| Massage heads rotate slowly. | Massage heads are too heavily loaded. | Release the massage heads. |
| Massage heads do not move. | Device not supplied with power. | Connect the device to the power supply. |

12. Warranty and service

Fitagon GmbH, Technoparkstrasse 1, CH-8406 Winterthur grants a warranty for this product under the following conditions and to the extent described below.

The following warranty conditions do not affect the seller's statutory warranty obligations arising from the purchase contract with the buyer. The guarantee also applies without prejudice to mandatory statutory liability provisions.

Fitagon guarantees the faultless functionality and completeness of this product. The warranty period is 2 years from the date of purchase of the new, unused product by the buyer.

This warranty only applies to products purchased by the buyer as a consumer and used exclusively for personal purposes within the scope of domestic use. Swiss law applies.

If, during the warranty period, this product proves to be incomplete or defective in function in accordance with the following provisions, Fitagon will provide a replacement delivery or repair free of charge in accordance with these warranty conditions.

If the buyer wishes to report a warranty claim, they should first contact Fitagon Customer Service:

Fitagon GmbH

Mail: info@fitagon.ch

Tel: +41 52 202 57 70

The buyer will then receive more detailed information on how to handle the warranty claim, e.g. where to send the product free of charge and which documents are required.

A claim under the guarantee can only be considered if the buyer

- a copy of the invoice/purchase receipt and
- the original product

Fitagon or an authorized Fitagon partner.

Expressly excluded from this guarantee are

- Wear and tear resulting from normal use or consumption of the product;
- Accessories supplied with this product that wear out or become damaged during normal use. are consumed (e.g. batteries, rechargeable batteries, cuffs, seals, electrodes, illuminants, attachments, inhaler accessories);
- Products that are used improperly and/or contrary to the provisions of the operating instructions have been used, cleaned, stored or maintained, as well as products that have been used, cleaned, stored or maintained by the buyer or not by Fitagon have been opened, repaired or converted;
- Damage that occurs during transport between the manufacturer and the customer or between the service center and the customer;
- Products that were purchased as 2nd choice items or as used items;
- consequential damage resulting from a defect in this product (however, in this case claims arising from product liability or other mandatory statutory liability provisions exist).

Repairs or a complete replacement do not extend the warranty period under any circumstances.

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